


Operating the Wave Software

To download the free Wave Software download it from this link:

http://www.synertronic.co.za/yahd/YahdProgInstaller_5_5.msi

Cut and paste this link
in your Windows
address Bar (Do not
paste it is Google
Search)

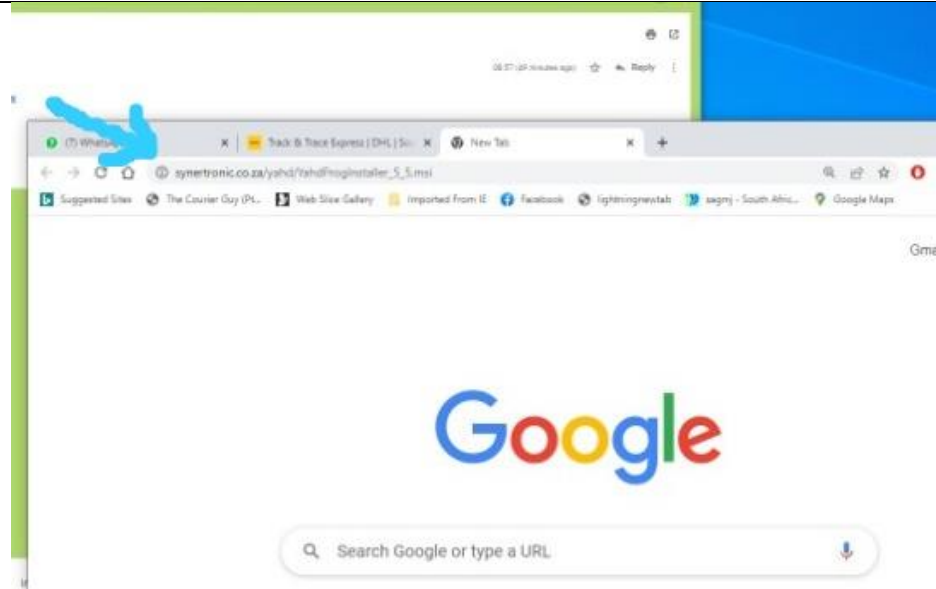

Click on “MORE INFO”
if you get this
message:






Click on “Run
Anyway”



Click on “Install” and follow the prompts to install the software.



Click on the icon and the INTERFACE screen will open.

1. Is the complete list
2. Is the default with trusted and tested treatments.
3. Use the “IMPORT LIST” to import updated lists that has been downloaded to your computer.

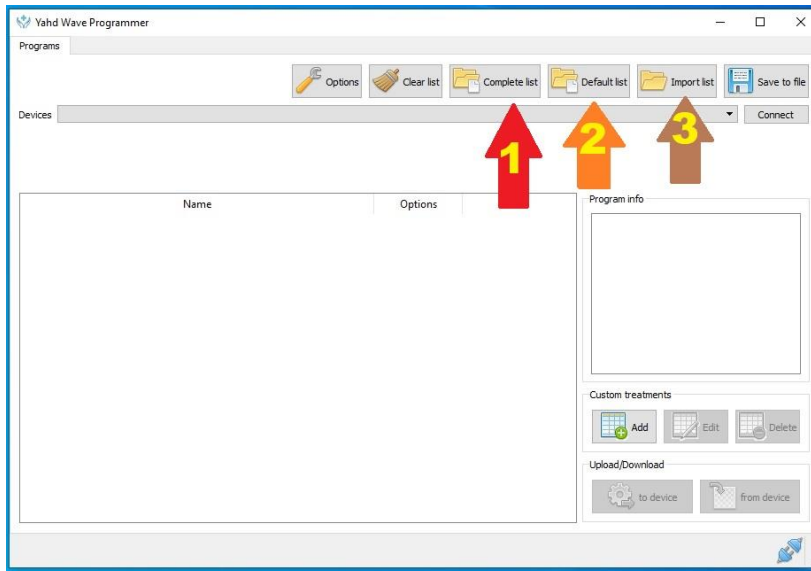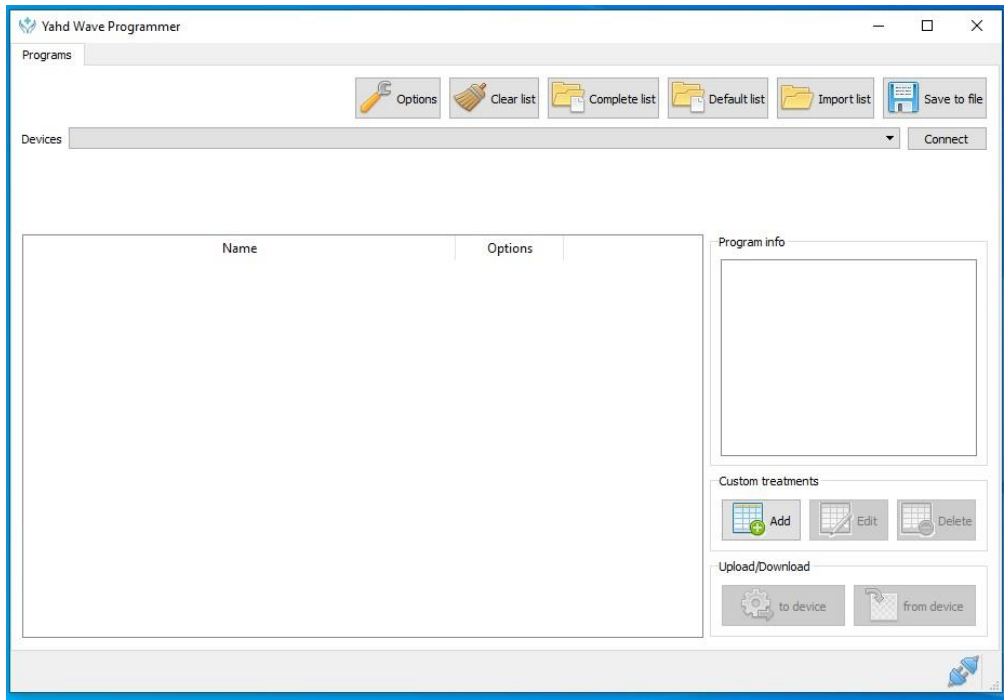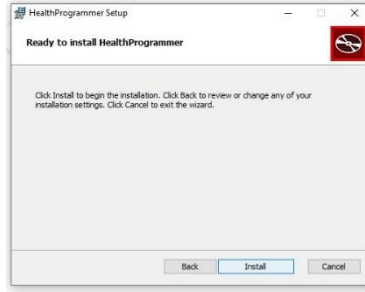



The INTERFACE screen with a uploaded selected list:

1. Click on "ADD" to add a Treatment or a SET OF TREATMENTS.
2. Select to add one treatment or to combine existing treatments on the list to a new SET.

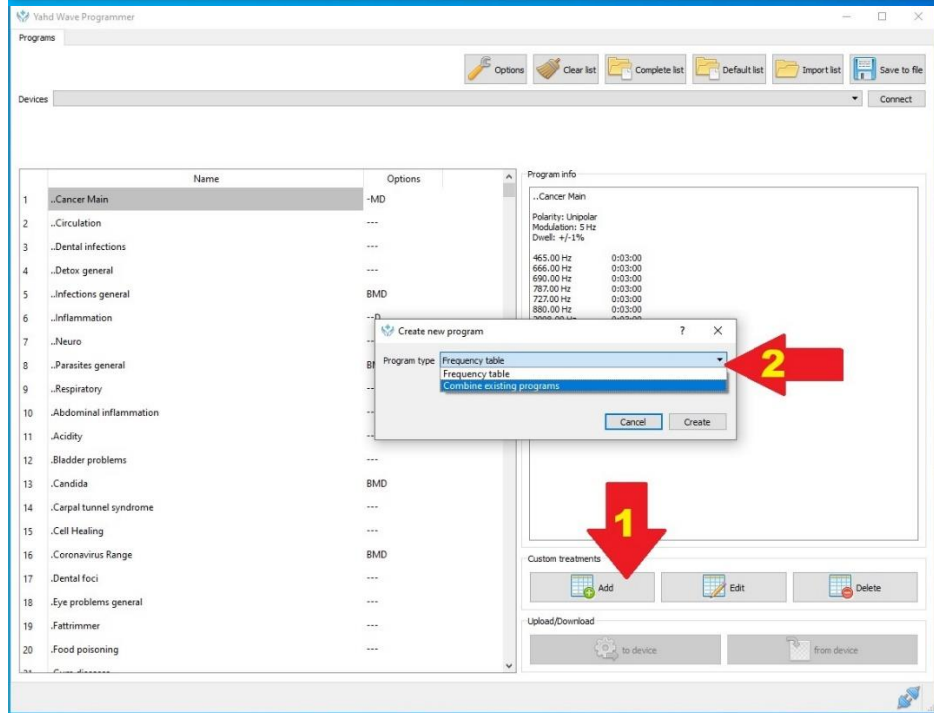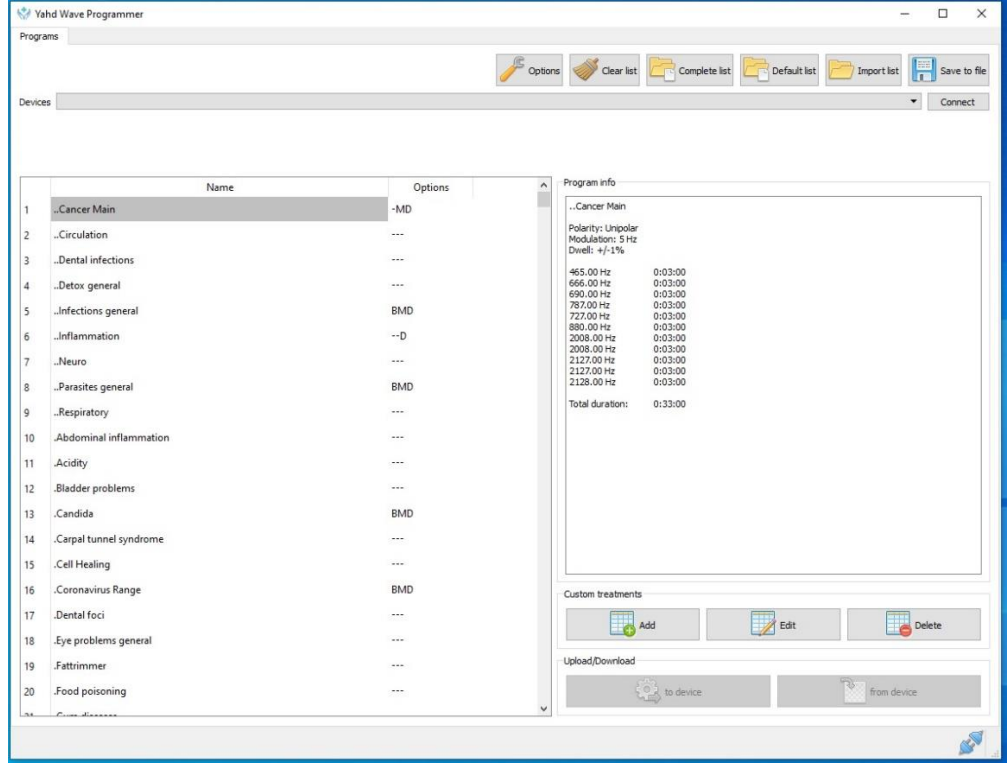


1. Give the new treatment or SET a new name.
2. Select the number of treatments that needs to be included in the SET.
3. Double click to activate the scroll to select a treatment from the list.

When task is completed click on "CREATE"

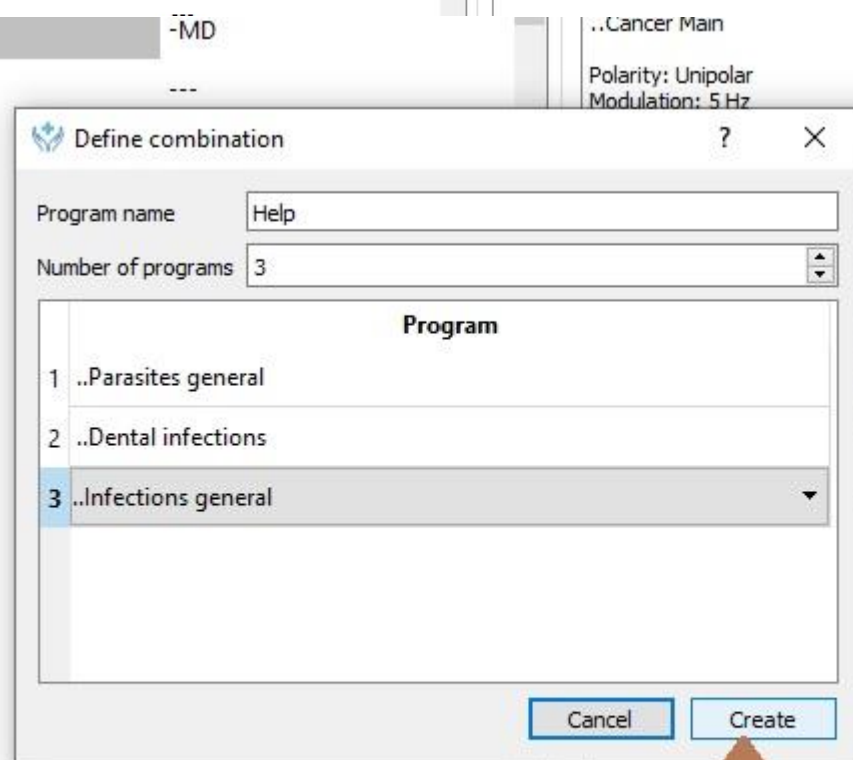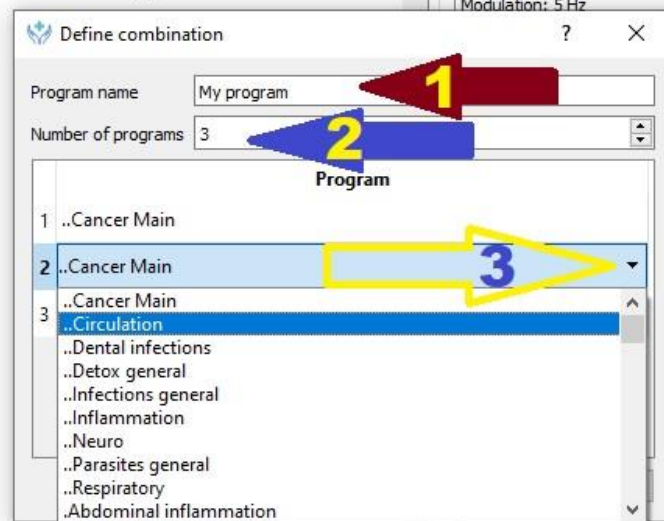

1. Your treatment or SET will be stored on the list that you have uploaded.
2. Save the adjusted treatment list to your computer.

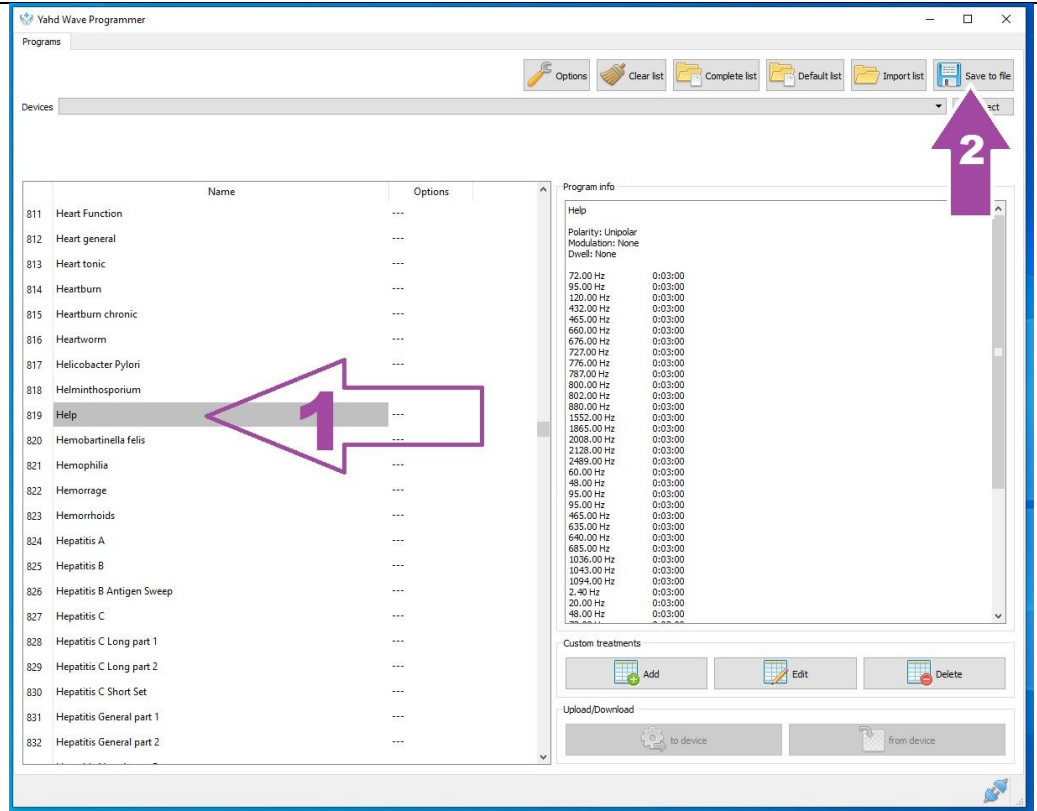

Connect your Wave Unit to your computer.



Click on "CONNECT"
on the INTERFACE
screen and your unit's
data will appear on
the top.
Click on "TO DEVICE"
to upload your own
updated list to your
unit.

